



Yes, you can.[®]

Invacare[®] Flo-tech[®] Contour

Designed for comfort and improved posture


The surface of the *Invacare[®] Flo-tech[®] Contour* cushion is contoured into sections. This unique design improves posture and increases the support surface area, effectively dispersing weight away from specific pressure points. The combination of improved posture and increased support surface area will contribute significantly to pressure ulcer prevention. This deep cushion is suitable for heavier clients (up to 30 stone), and those considered 'At Risk' (Waterlow score - up to 15)* of developing pressure ulcers.



3 year warranty

➤ Moulded foam cushion

Moulding allows **Flo-tech**® cushions to provide a 3D form, which cannot be emulated in cut foam. Moulded foam is highly durable, increasing the product's longevity. Theommel at the front of the cushion separates the thighs, encouraging good posture and allows air to circulate.

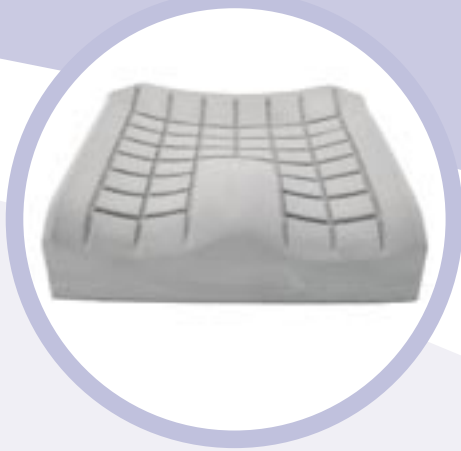

➤ Contoured seating surface

Offers improved stability and postural support. The castellations of the cushion allow slight movement which in turn reduces pressure. The contoured surface helps to prevent the client sliding forward, therefore reducing the potential shear and friction.

Options and Features



Water Resistant Cover
Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base.



Available in a wide range of sizes
The Invacare® **Flo-tech**® **Contour** is available in sizes: 14 x 13" to 24 x 22" (36 x 33 cm to 60 x 55 cm)*.


* See Invacare® **Flo-tech**® cushion size chart for details.

Technical data


Colours




Kingfisher Green PU (Standard) Black Smooth (Option) Black Wipedown (Option) Black Towelling (Option)

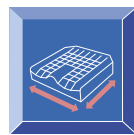

A wide range of sizes

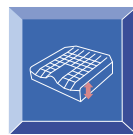


Height: 10.3 cm (4")

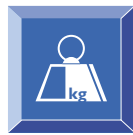

Weight: 1.38 kg



Max Client Weight: 191 kg (30 stone)

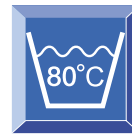

Cleaning: 80°C
Refer to instructions for use and relevant infection control guidelines.

Water Resistant
Vapour Permeable

Manufactured to European Standards: BS 7176 (Crib 5)

* Essential nursing care is pivotal in pressure ulcer prevention. These cushions will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited
Pencoed Technology Park · Pencoed · Bridgend CF35 5HZ
Tel: +44 1 656 776 222 · Fax: +44 1656 776 220
E-mail: uk@invacare.com · www.invacare.co.uk
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd
Unit 5 · Seatown Business Campus
Seatown Road · SWORDS · County Dublin
Tel. +353 1 8107084 · Fax +353 1 8107085
Email: Ireland@invacare.com · www.invacare.ie