



**Yes, you can.®**

## **Invacare® Flo-tech® Image**

### **Slimline contoured cushion with gel sac**

The *Invacare® Flo-tech® Image* cushion uses a combination of a slim, contoured foam with a large gel sac covering the sacral area. This increases even weight distribution and immersion of the Ischial Tuberosities. The design delivers a high degree of protection to clients who are considered to be at 'High Risk' of developing pressure ulcers (Waterlow score - 20)\*. This slim cushion format is ideal for both chair and wheelchair use.



**3 year warranty**

## ► Slimline contoured foam and gel combination cushion

Some clients prefer a slimline cushion, finding it more compatible with their daily living activities. The *Invacare® Flo-tech® Image* combines the benefit of a contoured foam upper surface with a large gel sac covering the ischial and sacral area, offering increased support and reducing the occurrence of peak pressures.



## ► Contoured seating surface with gel sac

Offers improved stability and postural support. The moulded foam castellations and gel combine effectively to reduce shear and friction forces on tissue, and help to prevent the client sliding forward in the seat.

### Options and Features



#### Water Resistant Cover

Cover is two-way stretch, water-resistant and vapour-permeable. Cover features an anti-slip base and carry handle.



#### Available in a wide range of sizes


The *Invacare® Flo-tech® Image* is available in sizes: 14 x 13" to 19 x 21" (36 x 33 to 48 x 53 cm). \*

\* See *Invacare® Flo-tech®* cushion size chart for details.


### Technical data

Warranty  
3 Years


#### Colours

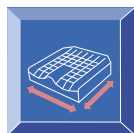

Wide range of sizes

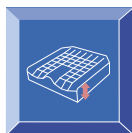

Height: 7.8 cm (3")

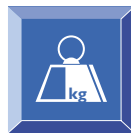

Weight: 1.82 kg



Max Client Weight: 140 kg (22 stone)



Cleaning: 80°C  
Refer to instructions for use and relevant infection control guidelines.

Cover:  
Water Resistant  
Vapour Permeable

Manufactured to European Standards:  
BS 7176 (Crib 5)



Mulberry PU (Standard)



Black Smooth (Option)



Black Wipedown (Option)



Black Towelling (Option)

\* Essential nursing care is pivotal in pressure ulcer prevention. These cushions will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

All information quoted is believed to be correct at time of print. *Invacare®* reserves the right to alter product specifications without prior consultation.



INVACARE® Limited  
Pencoed Technology Park · Pencoed · Bridgend CF35 5HZ  
Tel: +44 1 656 776 222 · Fax: +44 1656 776 220  
E-mail: uk@invacare.com · [www.invacare.co.uk](http://www.invacare.co.uk)  
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd  
Unit 5 · Seatown Business Campus  
Seatown Road · SWORDS · County Dublin  
Tel. +353 1 8107084 · Fax +353 1 8107085  
Email: Ireland@invacare.com · [www.invacare.ie](http://www.invacare.ie)