



**Yes, you can.<sup>®</sup>**

## **Invacare<sup>®</sup> Softform<sup>®</sup> Original**

### **The clinically proven replacement mattress**

The *Invacare<sup>®</sup> Softform<sup>®</sup> Original* is an established and clinically proven replacement mattress. Widely used across the UK healthcare sector, this high quality mattress offers unique design features to provide effective protection to those considered to be at 'High Risk' (Waterlow score - 20)\* of developing pressure ulcers.



### **Warranty**

**Foam: 4 Years**  
**Cover: 4 Years**

## ➤ Clinically proven pressure ulcer prevention


The *Invacare® Softform® Original* features dual contoured insert pads, constructed from 'keyhole cut' castellated foam cells, which move independently to maximise the body contact area and optimise pressure reduction.




## ➤ Flip and Rotate mattress

The robust, high density 'H' core ensures that the mattress does not bottom out and offers supportive side walls for transfers in and out of bed. In addition, the mattress can be flipped and rotated through a four-way turning cycle on a regular basis, effectively prolonging the mattress life.

## Features

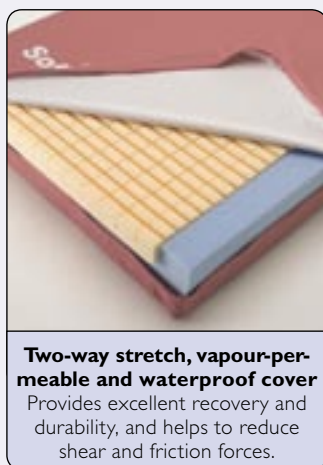

**Two-way stretch, vapour-permeable and waterproof cover**  
Provides excellent recovery and durability, and helps to reduce shear and friction forces.

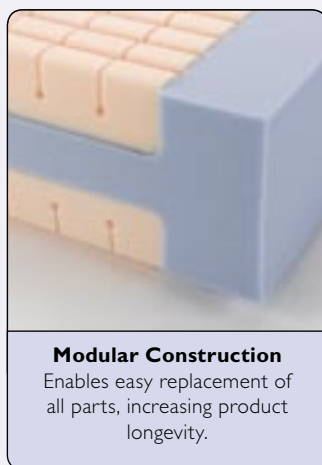


**Modular Construction**  
Enables easy replacement of all parts, increasing product longevity.




**Clinically proven**  
The *Invacare® Softform Original* has been proven in clinical trials and is an established choice with healthcare professionals.

## Technical data


**Warranty:**  
Cover: 4 Years  
Foam: 4 Years

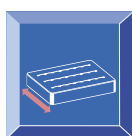

Std Width: 80, 85.5, 88, 90 cm

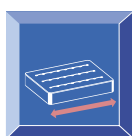

Std Length: 195 cm

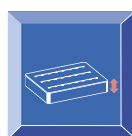

Height: 17.5 cm



Weight: 12.5 kg



Max. Client Weight: 26 stone (165 kg)



Cleaning: 80°C  
Refer to instructions for use and relevant infection control guidelines

## Colour



Pink (PU)

Additional sizes available on request

Manufactured to European Standards:  
BS 7177 Crib 5

\*Essential nursing care is pivotal in pressure ulcer prevention. This mattress will positively contribute to the outcome of a pressure ulcer prevention care plan. Education, clinical judgement and action based planning based on vulnerability are fundamental factors in the prevention of pressure ulcers. The Waterlow score is a formal method of assessing risk from pressure ulcer development, and should be used in conjunction with an informal assessment (informed nursing judgement). Informal assessment is considered to be of greater importance and clinical value.

All information quoted is believed to be correct at time of print. Invacare® reserves the right to alter product specifications without prior consultation.



INVACARE® Limited  
Pencoed Technology Park · Pencoed · Bridgend CF35 5HZ  
Tel: +44 1 656 776 222 · Fax: +44 1656 776 220  
E-mail: uk@invacare.com · [www.invacare.co.uk](http://www.invacare.co.uk)  
Sales Order E-mail: ordersuk@invacare.com

INVACARE® Ireland Ltd  
Unit 5 · Seatown Business Campus  
Seatown Road · SWORDS · County Dublin  
Tel. +353 1 8107084 · Fax +353 1 8107085  
Email: Ireland@invacare.com · [www.invacare.ie](http://www.invacare.ie)

