

QUANTUM
MOVING FORWARD

A Division of Pride Mobility Products® Ltd.



Athena

REAR-WHEEL DRIVE POWER CHAIR

Athena shown with Vogue seating, LED lighting package and optional 12" Lift and Tilt system

Unequaled Versatility


The Quantum® Athena rear-wheel drive power chair delivers aggressive, terrain-gripping torque and is available with high speed motors for a top speed of up to 10.46 km/h (6.5 mph). The innovative, compact design of the Quantum Athena makes it highly manoeuvrable in tight spaces without sacrificing traditional rear-wheel drive performance.



Specifications

WHEELS	Drive: 35.5 cm (14") pneumatic Casters: 23 cm (9") pneumatic (standard) Anti-tips: 4" solid
SUSPENSION	Sport Trac
MAXIMUM SPEED¹	Up to 10.46 km/h (6.5 mph)
BRAKES	Intelligent Braking (electronic, regenerative disc brakes)
GROUND CLEARANCE²	6.5 cm (2.625")
TURNING RADIUS²	60 cm (23.5")
OVERALL SIZE²	Length: 84 cm (33") Width: 65 cm (25.5")
SEATING	Synergy [®] Seat Adjustable Width: 25 cm - 51 cm (10"-20") Adjustable Depth: 25 cm - 51 cm (10"-20") TRU-Balance [®] Power Positioning System Adjustable Width: 35.5 cm - 51 cm (14"-20") Adjustable Depth: 35.5 cm - 51 cm (14"-20") Vogue Seating Adjustable Width: 41 cm - 51 cm (16"-20") Adjustable Depth: 41 cm - 51 cm (16"-20")
DRIVE TRAIN	Two-motor, rear-wheel drive
BATTERIES	2 x 12 Volt, deep cycle, 70 AH
BATTERY CHARGER	Off-board, 8 amp (standard)
AVAILABLE ELECTRONICS	70 amp, Pilot Plus (standard)
WEIGHT CAPACITY	136 kg (21 Stone) maximum (user weight)
BASE WEIGHT	65 kg (144 lbs.)
BATTERY WEIGHT³	70 AH: 24 kg (53.5 lbs)
SEAT WEIGHT⁴	20 kg (44 lbs.) (Synergy Seat)
WARRANTY	Lifetime limited warranty on frame 1-year limited warranty on electronics 18-month warranty on drive control system 2-year limited warranty on drive train

Shroud Colours

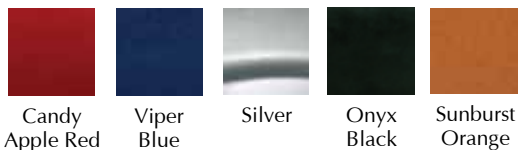

The information contained herein is correct at the time of publication; due to the manufacturers commitment to constant improvement and development, they reserve the right to alter specifications without prior notice.

Theoretical range is calculated under test conditions in accordance with European standards. Actual range in normal use will depend on many factors, including the condition of the vehicle and its batteries, the weight of the driver, correct tyre pressure, ambient temperature, and the gradient and surface of the road or pavement.



32 Wedgwood Road, Bicester,
Oxfordshire OX26 4UL
(Tel) 01869 324 600
(Fax) 01869 323 070
www.quantumrehab.co.uk

Features

- Easy operation results in a short learning curve
- Narrow 65 cm (25.5") width and compact 84 cm (33") length
- Outstanding tight-space manoeuvrability
- Sport-Trac Suspension
- Powerful motors for aggressive performance and terrain-gripping torque
- 3-point seat mount for ease of adjustment and addition of optional components
- Easy rear access to batteries
- Transit Securement Package allows user to conveniently secure Quantum[®] Athena for safe transport

Options

- Synergy[®] seating system
- TRU-Balance[®] power positioning system
- TRU-Comfort[™] seating system
- Lap belt
- Flat-free tyres
- Cup holder
- Vogue seating
- 70° swingaway footrests
- Swing-away joystick mount
- 80° heavy-duty leg rests
- Elevating leg rests
- Heel loops
- Oxygen tank holder
- Cane/crutch holder
- Rear basket

FOOTNOTES:

- 1) Speed and range vary with user weight, terrain type, battery charge, battery condition and tyre pressure
- 2) Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3 %.
- 3) Battery weight may vary based on manufacturer and battery size
- 4) Seat weight is based on an 46 cm (18") Synergy seat and will vary by seat size, specification and options chosen.

